





Sleeve Installation Guide

Hip Pro Smart Wrap, Right


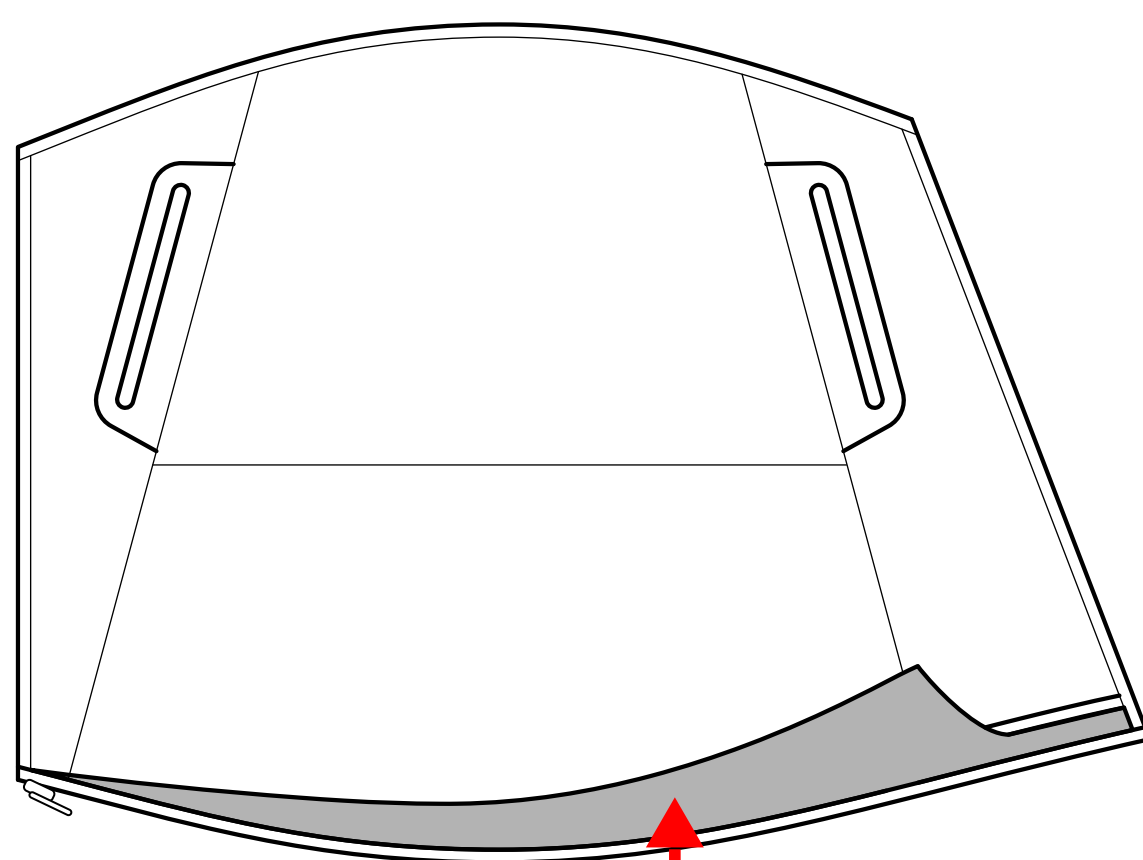


Right Hip Sleeve

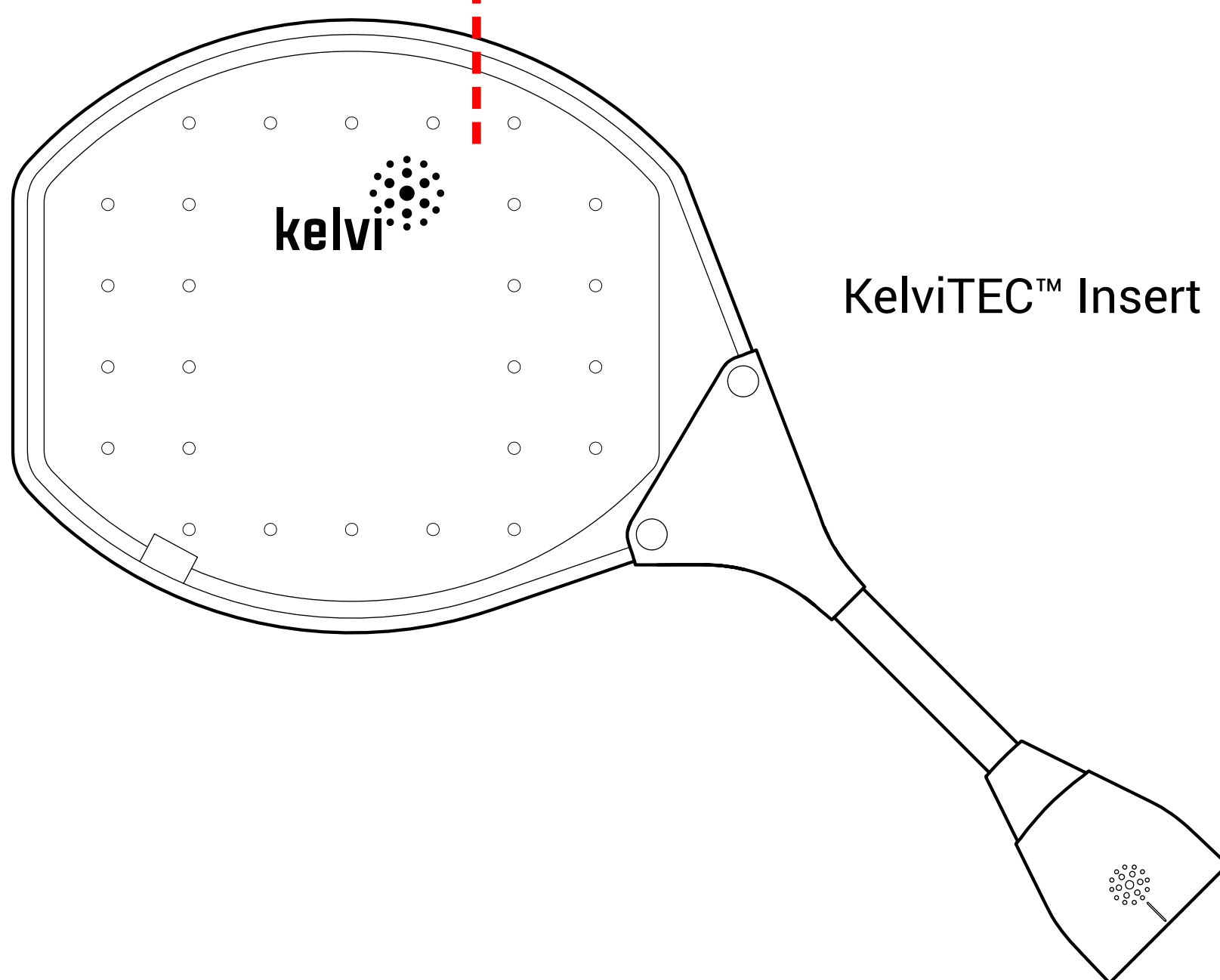

KelviTEC™ Insert

1. Place the Hip Sleeve on a flat surface and open the zipper.
2. Insert the KelviTEC Insert™ into the sleeve with the metal tiles facing down, away from the Velcro side.
3. Ensure the KelviTEC Insert™ is fully inserted and close the zipper.

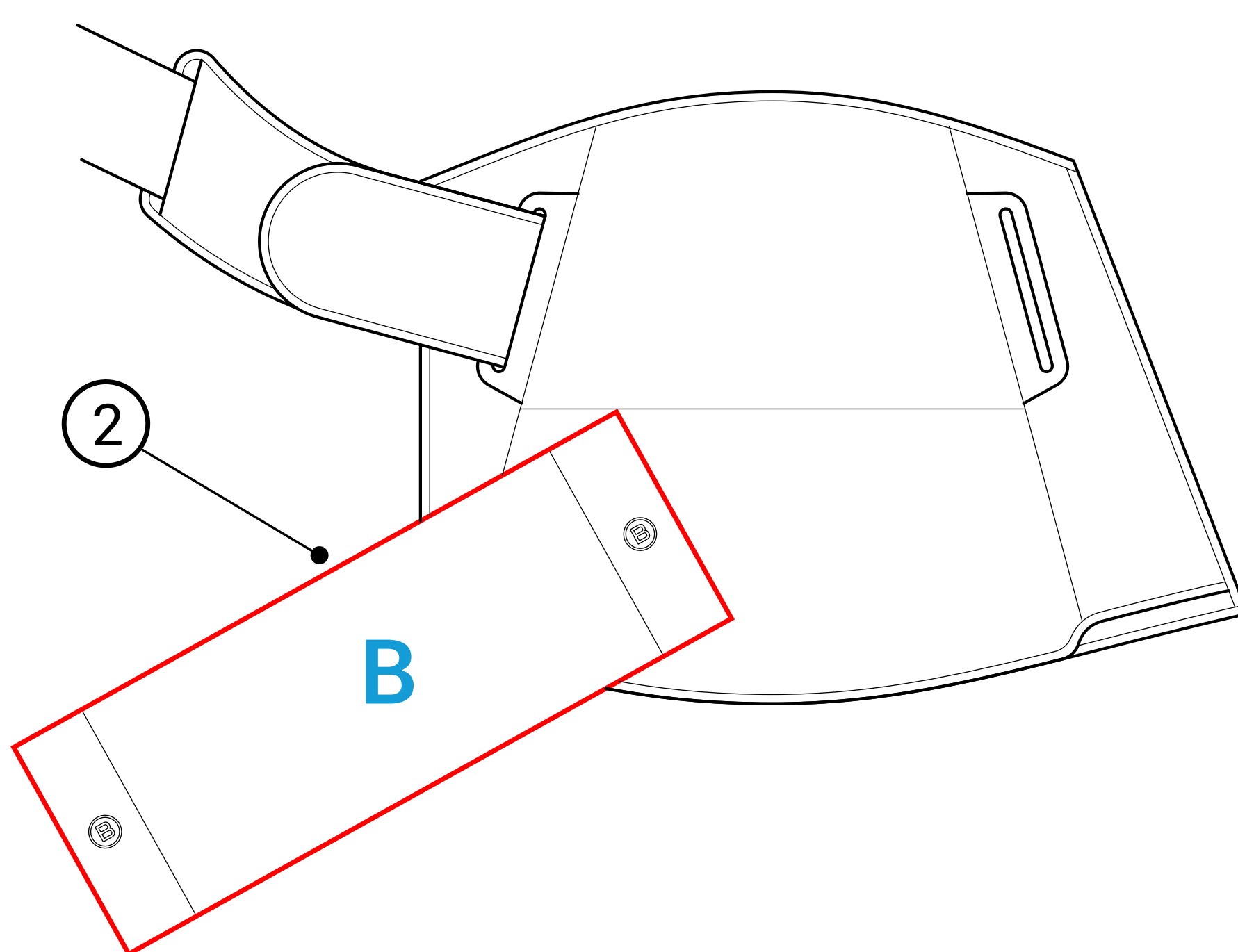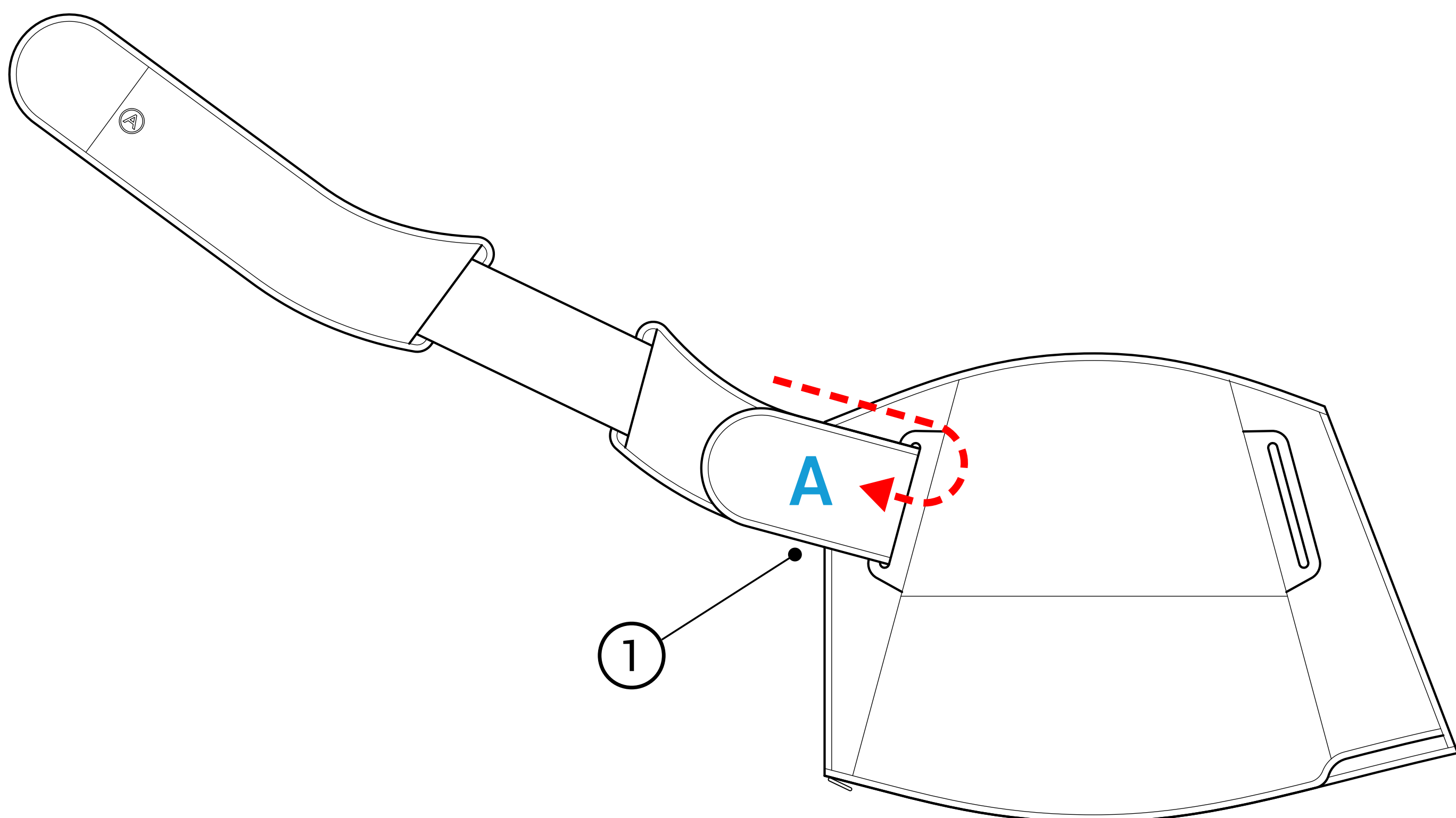

Prepare to Wear

1. Attach the waist strap **A** to the left D-ring, leaving the opposite side unattached.
2. Attach the thigh strap **B** as shown.