





Sleeve Installation Guide

Shoulder Pro Smart Wrap

Smart Wrap™ Assembly

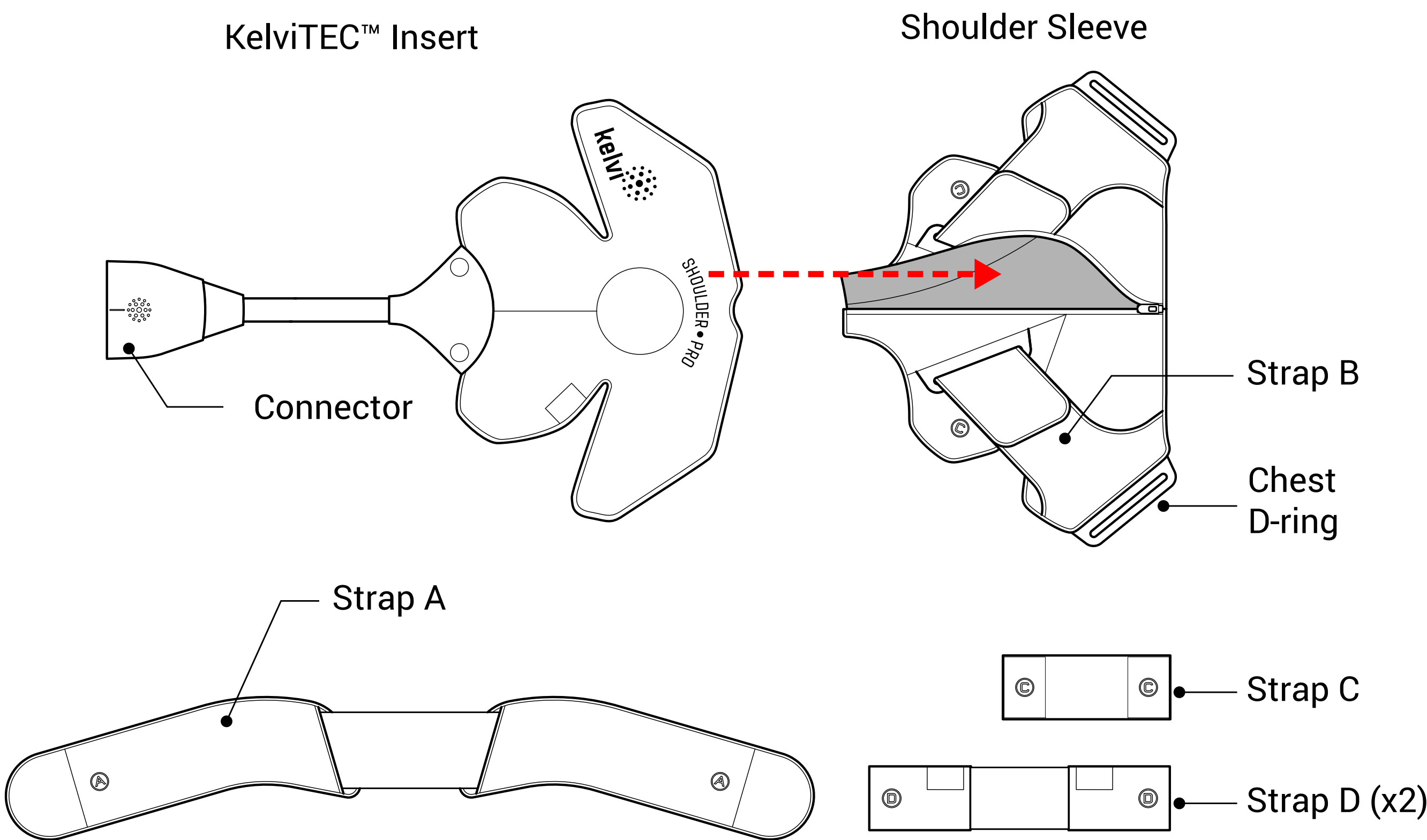

1. Place the Shoulder Sleeve on a flat surface and open the central zipper.
2. Insert the KelviTEC Insert™ into the sleeve with the metal tiles facing down, away from the zipper.
3. Ensure the KelviTEC Insert™ is fully inserted and close the zipper.

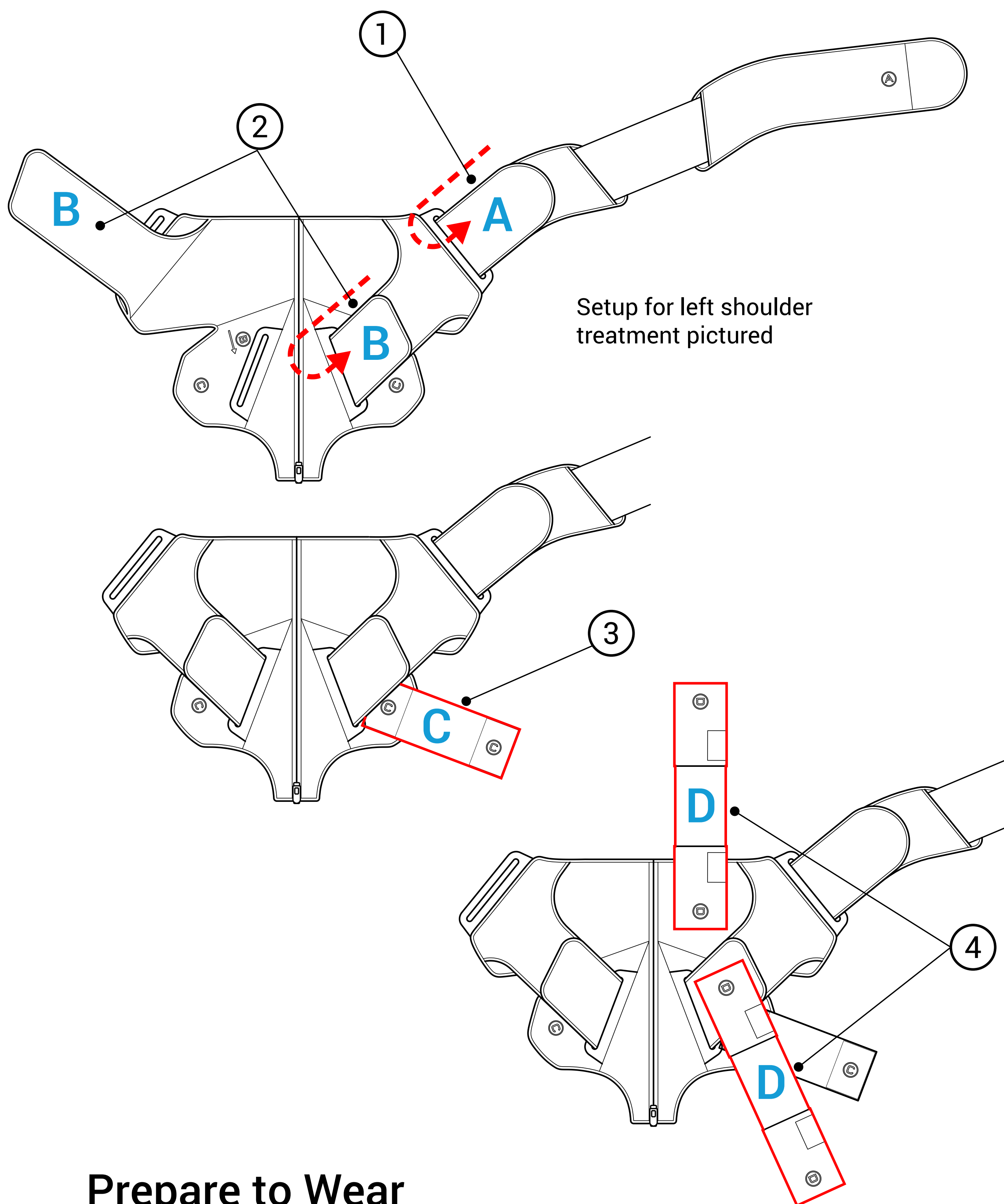

Prepare to Wear

1. Attach the chest strap **A** to the D-ring, leaving the opposite side unattached. For left shoulder treatment, attach the strap to the right D-ring. For right shoulder treatment, attach the strap to the left D-ring.
2. Insert both **B** straps through the indicated D-rings.
3. Attach the bicep strap **C** as shown
4. Attach the 2 auxiliary straps **D** as shown.